
beper

IT	CUSCINO MASSAGGIANTE	pag. 2
EN	MASSAGER PILLOW	pag. 5
FR	COUSSIN MASSANT	pag. 8
DE	MASSAGE-KISSEN	pag. 11
ES	COJÍN MASAJEADOR	pag. 14
GR	ΜΑΞΙΛΑΡΙΑ ΜΑΣΑΖ	pag. 17
RO	PERNĂ CU MASAJ	pag. 20
CZ	MASÁŽNÍ POLŠTÁŘ	pag. 23
NL	MASSAGEKUSSEN	pag. 26
LV	MASĀŽAS SPILVENS	pag. 30

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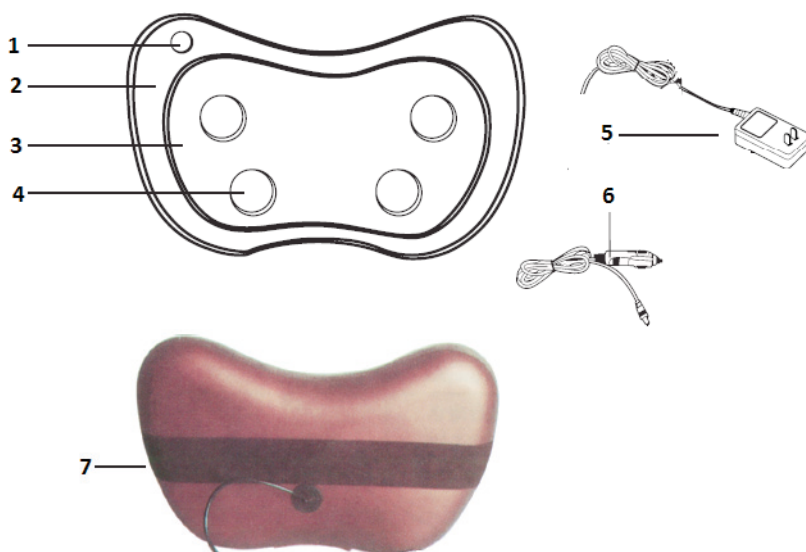assistenza@beper.com

Precautions

Read the following instructions carefully, since they give you useful safety information about installation, use and maintenance and help you avoid mishaps and possible accidents.

Remove the packing and make sure that the appliance is intact, with special attention to the supply cable.

The elements of the packing (plastic bags, polystyrene, etc) do not have to be accessible to the children, in order to preserve them from danger sources. It is recommended not to throw these packing material in the domestic waste, but to deliver them to the appropriate station of collection or destruction refusals, asking eventually information to the manager of the Sanitation Department of your own municipality.

If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.

Before connecting the appliance make sure that the mains voltage corresponds to the data specified on the label.

The appliance can be damaged if the voltage is not correct.

If the plug does not suit the socket, let it be replaced with another suitable plug by qualified personnel only, which shall also ensure that the section of the socket cables is suitable for your appliance power absorption.

Generally, the use adaptors multiple plugs and/or extensions is not recommended; should they be necessary, please use only simple or multiple adaptors and extensions in compliance with the safety rules in force, making sure that the capacity absorption limit marked on the simple adaptors and the extensions and the maximum power limit marked on the multiple adaptor are not exceeded.

Do not use the appliance if the cable or the plug shows any damages. Do not use the appliance after a possible irregular working. In this case, switch it off and do not break it open. Refer to an authorized Beper after-sale service for repair and require original spare parts only. Failing to follow these instructions may impair the safety of your appliance.

This appliance is to be destined to its special purpose only. Any other use is to be considered as improper and thus dangerous. The manufacturer cannot be responsible for any damages deriving from improper, wrong or incautious use.

This appliance is for household use only.

To protect against electric shock, do not immerse cord, plug or any parts of the appliance into water or other fluids.

Do not use the appliance near explosives, highly flammable materials, gases or burning flames.

Do not allow children, disabled people or every people without experience or technical ability to use the appliance if they are not properly watched over. Close supervision is necessary when any appliance used by or near children.

Components Fig.A

1. ON/OFF button
2. Synthetic leather cover
3. Massaging heads
4. Massaging wheels
5. Home adapter
6. Car adapter
7. Elastic band

Use

The massager pillow can be used on different parts of the body. It is suitable for neck, back, stomach, legs, arms, ankles and it helps against muscle contracture, besides promoting a better blood circulation.

On the back of the massager pillow there is an elastic band which permits to settle it. Place the massager pillow on the desired area and press the ON/OFF button to start it.

The massager pillow can be used at home or in the car, thanks to the included adapters.

After using the device, press again the ON/OFF button to stop it.

Attention

Don't use the massager pillow for longer than 15 minutes. After this time, the massager will shut off automatically.

Don't use the massager pillow on the stomach area directly after meals, as it could interfere with the digestion.

The massager pillow is not a medical device: in case of muscle contracture or other problems always talk with your doctor first.

Don't immerse the product in the water: only use a damp cloth to clean it. Don't use chemical detergents.

Technical Data

Power: 20 W

Power supply: 220V – 12V 2A

For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment.

The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

GUARANTEE CONDITIONS

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product

Contact your distributor in your country or after sales department Beper.

E-mail assistenza@beper.com which will forward your enquiries to your distributor.

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⚠ Attention

Don't use the massager pillow for longer than 15 minutes. After this time, the massager will shut off automatically. Don't use the massager pillow on the stomach area directly after meals, as it could interfere with the digestion.

The massager pillow is not a medical device: in case of muscle contracture or other problems always talk with your doctor first. Don't immerse the product in the water: only use a damp cloth to clean it. Don't use chemical detergents.

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