
beper

IT	SPAZZOLINO DA DENTI ELETTRICO CON CLESSIDRA	pag. 2
EN	ELECTRIC TOOTHBRUSH WITH HOURGLASS	pag. 4
FR	BROSSE À DENTS ÉLECTRIQUE AVEC SABLIER	pag. 6
DE	ELEKTRISCHE ZAHNBÜRSTE MIT SANDUHR	pag. 9
ES	CEPILLO DE DIENTES ELÉCTRICO CON RELOJ DE ARENA	pag. 11
GR	ΗΛΕΚΤΡΙΚΗ ΟΔΟΝΤΟΒΟΥΡΤΣΑ ΜΕ ΚΛΕΨΥΔΡΑ	pag. 14
RO	PERIUȚĂ DE DINȚI ELECTRICĂ CU CLEPSIDRĂ	pag. 17
CZ	ELEKTRICKÝ ZUBNÍ KARTÁČEK S PŘESÝPACÍMI HODINAMI	pag. 19
NL	ELEKTRISCHE TANDENBORSTEL MET ZANDLOPER	pag. 22
LV	ELEKTRISKĀ ZOBU BIRSTE AR SMILŠU PULKSTENI	pag. 24

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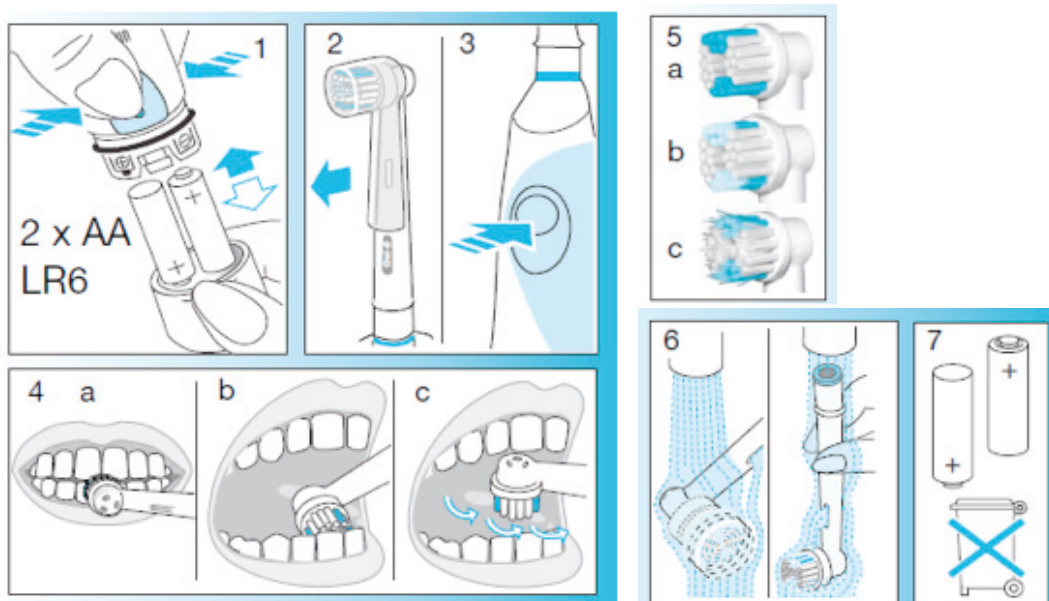assistenza@beper.com

General safety advice

Read the operating instructions carefully before first use of the appliance.

Read the following instructions carefully, since they give you useful safety information about installation, use and maintenance and help you avoid mishaps and possible accidents.

Remove the packing and make sure that the appliance is intact, with special attention to the supply cable. The elements of the packing (plastic bags, polystyrene, etc) do not have to be accessible to the children, in order to preserve them from danger sources.

It is recommended not to throw these packing material in the domestic waste, but to deliver them to the appropriate station of collection or destruction refusals, asking eventually information to the manager of the sanitation department of your own municipality.

This appliance is to be destined to its special purpose only.

Any other use is to be considered as improper and thus dangerous. The manufacturer cannot be responsible for any damages deriving from improper, wrong or incautious use.

To protect against electric shock, do not immerse any parts of the appliance into water or other fluids.

Do not expose the appliance to atmospheric agents such as rain, wind, snow, hail.

Do not use the appliance near explosives, highly flammable materials, gases, burning flames, heaters.

Do not allow children, disabled people or every people without experience or technical ability to use the appliance if they are not properly watched over. Close supervision is necessary when any appliance used by or near children.

Make sure that children do not play with this appliance.

If you decide not to use the appliance any longer, make it unserviceable.

Also make any dangerous parts unserviceable, especially for children who could play with the appliance.

Instruction for use

Insert batteries as shown. Uses two 1.5 V type AA (LR 6).

Wet brushhead and apply toothpaste.

Place brushhead on teeth and switch on.

Brush for at least 3 minutes (recommended by dentists). Guide the brushhead from tooth to tooth, first the outsides (a), then the insides (b), finally the chewing surfaces (c).

During the first days of use, gums may bleed slightly. Should bleeding persist after 2 weeks, contact your dentist.

Replace brushhead when indicated (a = new, b/c = replace)

Clean as shown under running water, then dry.

Dispose of discharged batteries according to national or local regulations.

Adult supervision recommended when used by children.

Technical data

Use 2xAA batteries (not included)

For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment.

The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guaranteee.

The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

GUARANTEE CONDITIONS

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product

Contact your distributor in your country or after sales department beper.

E-mail assistenza@beper.com which will forward your enquiries to your distributor.

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