Cod.: P302VIS003

Fig.A

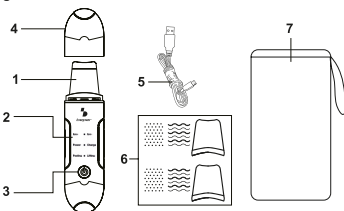
Fig.B

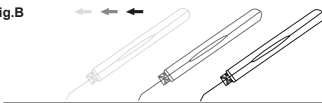
Fig.C


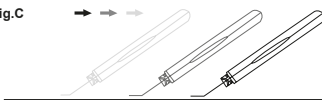
Fig.D

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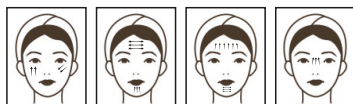EN - Skin scrubber - ultrasonic face cleaning device

Description of product Fig.A

1. Steel shovel
2. Operating mode LED indicator lights
3. ON/OFF key
4. Steel shovel protective cap
5. USB cable
6. N.2 silicone accessories for massage/cleansing
7. Storing bag

Keep the steel shovel at 45° during use as shown in fig. B and C; it is advisable not to insist on the same area of the face for long, but to gently slide the shovel on all the face area, as shown in fig. D. In case of skin irritations/inflammations, do not use the device and wait until the symptoms disappear.

For best results, it is advisable to use the device for maximum 2-3 times a week, do not use it daily and/or multiple times a day.

Use Fig. B-D

Before first use, charge the battery completely, by inserting the provided USB cable inside the connector found at the base of the device. Connect the device to a wall socket (via a USB adapter not included), or to other charging devices (e.g. laptop, powerbank, PC, etc.)

Once connected to the charging source, all the LED indicators on the device will start blinking, and will become fixed once battery is fully charged.

When the battery is dead, the device will buzz and will automatically shut off; charge the battery to use the device again. To switch the device on, hold the ON/OFF key for a few seconds, the device will buzz and will turn on in "Ultrasonic cleaning" mode (Peeling • Lifting) – the corresponding LED indicator will light up.

To change mode: press the ON/OFF key once, to activate "EMS massage" mode (Power • Charge), the LED indicator will flash during operation; press the ON/OFF key twice to activate "Ion+" mode (Ion+), the corresponding LED indicator will light up; press the ON/OFF key three times to activate "Ion-" mode (Ion-), the corresponding LED indicator will light up; press the ON/OFF key four times to go back to "Ultrasonic cleaning" mode (Peeling • Lifting). To switch the device off, hold the ON/OFF key for a few seconds, the device will buzz and all the LED indicators will turn off.

Operation modes

Choose the operation mode based on the treatment you wish to carry out, and follow the suggested paths in fig. D based on the area to be treated.

Ultrasonic cleaning (Peeling • Lifting)

This operation mode allows a deep face cleaning, removing impurities and contributing to make skin smoother.

Before starting, slightly wet your face with warm water and apply a face cleansing solution, and rub gently to form a light foam; for more effective cleaning, place the cleansing or massaging silicone accessories on the shovel. When your face (or area to be treated) is covered in foam, remove it with the device making sure to move the shovel as shown in fig. B (if present, remove the silicone accessory before proceeding).

Proceed for about 3 minutes, gently moving the device all over your face. Finally rinse your face with warm water.

Ion+

This mode is ideal for a deep cleansing of face pores. Apply a cleansing lotion/toner on the skin, and massage making sure to move the shovel as shown in fig. B; after about 3 minutes, rinse your face with warm water.

Important: do not use the silicone accessories for cleansing or massaging, to avoid interfering with the ion export process and compromising the efficacy of treatment.

Ion-

This mode enhances the absorption of nutrients by the skin, improving skin hydration.

After properly cleansing your face, apply a moisturizing cream and use the device making sure to move the shovel as shown in fig. C, for about 3 minutes.

Important: do not use the silicone accessories for cleansing or massaging, to avoid interfering with the ion export process and compromising the efficacy of treatment.

EMS massage (Power • Charge)

This mode allows to carry out a gentle massage to enhance skin tone. After properly cleansing your face, leave the skin slightly wet and use the device making sure to direct the shovel as shown in fig. C for about 3 minutes.

Cleaning and maintenance

Keep this device away from water and excessive humidity. Make sure no part of the device (except the shovel and the silicone accessories) comes into contact with water or other liquids. Before cleaning, make sure the device is switched off. Use a soft cloth slightly dampened with water to remove any stain and/or residues of lotion from the shovel and the device. After each use, thoroughly wash the shovel and the silicone accessories with mild neutral soap or other similar detergents. Do not rinse the device and the shovel directly under tap water, to avoid damaging the device.

If the device is not used for long periods of time (approx. 3 months or more), it is advisable to fully charge the battery to avoid damaging it.

Technical data

3,7V - 0.65A 1.5W
Lithium battery 500 mAh
Input: 5V DC USB

For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream.

Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment. The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

Contact your distributor in your country or after sales department beper.

E-mail assistenza@beper.com which will forward your enquiries to your distributor.

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