
MASCHERA VISO A LUCI LED	IT
LED LIGHT FACE MASK	EN
LED-LICHT-GESICHTSMASKE	DE
MASQUE FACIAL À LUMIÈRE LED	FR
MÁSCARA FACIAL CON LUZ LED	ES
LED MASKA ZA LICE	HR
LED MASKA ZA OBRAZ	SL
ΜΑΣΚΑ ΠΡΟΣΩΠΟΥ ΜΕ LED ΦΩΣ	GR
MASCA FACIALĂ CU LUMINĂ LED	RO
LED SVĚTELNÁ MASKA NA OBLIČEJ	CZ
LED SEJAS MASKA	LV
LED VEIDO KAUKĖ	LT
LED-VALGUSEGA NÄOMASK	EE
LED-GEZICHTSMASKER	NL
LED MASKA ZA LIČE	BG
LED MASKA ZA LICE	RS
LED-FÉNYES ARCPRAKOLÁS	HU

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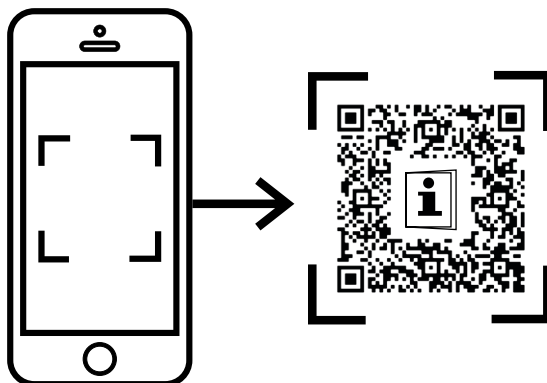## **General warnings**

### **Read the following instructions before using the appliance.**

Some elementary precautions must be followed before and during use of the appliance.

After removing the packaging make sure of the integrity of the appliance. In case of doubt do not use the appliance and contact professionally qualified personnel. Packaging elements (plastic bags, polystyrene, etc.) should not be left within the reach of children, as they are potential sources of danger.

Always make sure that the mains voltage is the same as that indicated on the technical data label and that the system is compatible with the power of the appliance.

Never pull the cord to disconnect it from the electrical outlet.

Make sure the cable is not in contact with hot or sharp surfaces.

Do not plug the appliance into the electrical outlet for charging the battery if the cable is damaged.

If the charging cable of the appliance is damaged, it must be replaced by the Manufacturer or his Technical Assistance Service or in any case by a person with a similar qualification, in order to prevent any risk.

In general, the use of adapters, multiple sockets and/or extension cords is not recommended. If their use becomes essential, it is necessary to use only adapters and extensions that comply with current safety regulations.

This appliance shall be intended only for the use for which it was expressly designed. Any other use is to be considered improper and therefore dangerous, as well as to cause the warranty to lapse. The manufacturer cannot be considered liable for any damage resulting from improper, erroneous and unreasonable uses.

To avoid dangerous overheating, it is recommended to unwind the power cord along its entire length and unplug it from the power supply network when the appliance is not in use.

Before carrying out any cleaning or maintenance operation, disconnect the appliance from the power supply network.

Never immerse the appliance in water or other liquids.

Do not keep the appliance near heat sources (e.g. radiator).

This appliance should not be used by people (including children) with reduced physical, sensory or mental abilities; by people who lack experience and knowledge of the appliance, unless they are carefully

supervised or well educated regarding the use of the appliance itself by a person responsible for their safety.

Make sure children do not play with braces.

When it is decided to dispose of this appliance as waste, it is recommended that it be rendered inoperative. It is also recommended to render harmless those parts of the appliance which are likely to constitute a hazard. This appliance should not be used by children. Keep the appliance and its cable out of the reach of children.

**Children must not play with the appliance. Keep the appliance and its cable away from children under 8 years of age.**

**Never immerse the motor body in water or other liquids. Never touch moving parts of the appliance.**

**⚠ Attention: Battery life and performance may be affected by several factors, including time of use, environmental conditions, and charge and discharge cycles. With the passage of time and prolonged use of the product, there may be a natural decrease in battery performance.**

### **Warnings for use**

Do not immerse the device, USB cable or remote control in water. Keep away from liquids and/or heat sources.

Do not remove components, screws, etc. and do not attempt to disassemble or repair the appliance.

Do not use the device for more than 15 minutes of treatment (maximum timer duration available).

Do not use the product when charging.

**⚠ Please note: This device is designed for use on adults over the age of 18. The product cannot be used by pregnant women, children, people with epilepsy, photosensitive people, or people using antihistamines and photosensitive drugs.**

This appliance may not be used by people with reduced physical, sensory or cognitive abilities, or with little experience and knowledge, without supervision or without having been given indications and warnings as to the correct use of the apparatus and the possible risks involved.

Do not use the device on skin conditions, skin ulcers, traumatic or post-operative wounds; before using it, consult your doctor if you have any physical conditions or conditions that require special attention (e.g., wounds, muscle trauma, etc.).

This product is not a medical device, and is not intended for the diagnosis, treatment or prevention of diseases and/or pathologies.

Do not look directly into the lights as this may cause temporary discomfort. Be sure to place the device on your face so that you can see through the eye holes and that the device does not cover your eyes.

Do not use if the skin is red or swollen due to previous treatment.

Discontinue use if the device overheats or becomes uncomfortably hot.

If you have any concerns or concerns about the skin you are treating or any medical questions regarding the use of this device, stop using it and consult your doctor.

#### **Product Description Fig. A**

1. Flexible silicone face mask
2. USB Type C Power Cable
3. Eye-darkening inserts
4. Detachable command
5. Adjustable bands with velcro
6. Type-C charging cable
7. LED light points
8. Timer selector
9. Light Mode Selector
10. ON/OFF button
11. USB charging socket

#### **Charging**

Before first use, you need to fully recharge the appliance.

Connect the USB cable to the controller's charging socket and then to a USB charging source.

While charging, the LED lights will start flashing. Once charging is complete, the LED lights will remain fixed and on. At this point you can disconnect the USB cable from the controller.

A full charge lasts approximately 3 hours, which varies depending on the power of the USB charging source. A full charge allows up to approximately 1 hour of use with cycles of approximately 15 minutes each.

The appliance switches off autonomously at the end of an operating cycle of the selected duration. You can change modes during treatment and the timer continues with its preset cycle.

#### **Light sensitivity test**

Before first use, perform a light sensitivity test.

Turn on the device and place it on one arm, allowing it to radiate for 10 minutes. If redness, itching or other symptoms appear during treatment, or within two hours of its conclusion, this may indicate skin photosensitivity. It is recommended to discontinue use and consult a dermatologist before proceeding with further applications.

Before use, clean your face and keep it completely dry (without makeup), as any cosmetic containing photosensitizing ingredients may reduce or affect the effectiveness of the treatment.

## Use

The product is designed exclusively to restore the skin's radiance and elasticity. It also helps improve skin metabolism, promote anti-aging action, and reduce inflammation, acne, and wrinkles.

Remove the product from the package and check that all parts are included and intact.

Insert the adjustable bands at the side ends of the mask.

The device is safe for use even without eye-darkening inserts. If the LED lights turn out to be too bright, add the inserts at the eye openings for a more comfortable experience. On the back of the blackout inserts are the letters "R" and "L", which indicate the insert for the right eye (Right) and the one for the left eye (Left) respectively. This marked part must face the inside of the LED mask, corresponding to the area with the light points. Be sure to position them correctly before use.

Bring the mask closer to your face and join the ends of the bands behind your head, adjusting them as needed.

Connect the USB power cable of the LED mask to the controller.

Press and hold the ON/OFF button for a few seconds. The device turns on and automatically sets to red light.

When the device is activated, press the timer selection button consecutively to set a usage cycle of 5, 10 or 15 minutes.

**⚠ Caution: Do not use the device for more than 15 minutes of treatment.**

You can select the desired light mode using the LED light mode selector.

By pressing the LED light mode button consecutively, the lights activate and allow the following different modes:

Mode	Wavelength	Indicated for
Red Light	650 nm	stimulates collagen regeneration, reduces wrinkles and improves skin elasticity
Orange Light	605 nm	promotes blood circulation and accelerates metabolism to repair damaged skin and reduce pigmentation
White Light	400-550 nm	accelerates skin renewal, helping to lighten spots and promoting skin relaxation
Purple Light	430 nm	dual-frequency combination of red and blue light, acting as an anti-acne effect and skin repair
Green Light	530 nm	illuminates and rebalances water and sebum levels in the skin, with a stress and tension relief effect
Blue Light	470 nm	promotes blood circulation and reduces acne inflammation
Yellow Light	590 nm	evens out skin tone for a more even, clean and bright complexion

At the end of the preset cycle, the device automatically turns off. If you want to stop treatment early, simply hold down the ON/OFF button for a few seconds to turn off the device.

Remove from the face and place in a dry, clean place. Proceed with your usual facial care routine.

For best results, regular use of the appliance is recommended. The intensity and speed of the results may vary depending on personal conditions.

After use, your skin may have a slight redness, a normal skin reaction that will return to normal within 1-2 hours.

**⚠ Caution: Stop using immediately if you experience burning or other abnormal discomfort during or after treatment.**

## Cleaning and maintenance

Clean the device after each use. Clean the outer surface with a soft, slightly umi-do cloth. Do not use detergents or abrasive solutions, these may damage the product coating.

All parts must be completely dry before storing the mask.

Store the product in a dry, clean place away from heat and direct sunlight. Store the product in the original packaging for maximum protection.

#### Technical data

Rechargeable Li-ion Battery 3.7V-800mAh

Input: 5V USB

**With a view to continuous improvement, Beper reserves the right to make changes and improvements to the product in question without prior notice.**

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment.

The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

#### GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

**The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance**

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

#### GUARANTEE CONDITIONS

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product

Contact your distributor in your country or after sales department Beper.

**E-mail assistenza@beper.com** which will forward your enquiries to your distributor.

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