
# beper

IT	BILANCIA PESAPERSONE IMPEDENZIOMETRICA CON APP	pag. 2
EN	BIOELECTRICAL IMPEDANCE BODY SCALE WITH APP	pag. 6
FR	BALANCE DE PESAGE IMPEDENZIOMETRIQUE AVEC APP	pag. 10
DE	KÖRPERANALYSEPERSONENWAAGE MIT APP	pag. 14
ES	BÁSCULA DE BAÑO DE IMPEDANCIA CON APLICACIÓN	pag. 18
GR	ΖΥΓΑΡΙΑ ΣΩΜΑΤΟΣ ΒΙΟΗΛΕΚΤΡΙΚΗΣ ΑΝΤΙΣΤΑΣΗΣ ΜΕ APP	pag. 22
RO	CÂNTAR CORPORAL CU IMPEDANȚĂ BIOELECTRICĂ CU APLICAȚIE	pag. 26
CZ	VÁHA - TĚLESNÝ ANALYZÉR S BT A APP	pag. 30
NL	BIO-ELEKTRISCHE IMPEDANTIE LICHAAMSWEEGSCHAAL MET APP	pag. 34
LV	BIOELEKTRISKĀS PRETESTĪBAS ĶERMEŅA SVARI AR APP	pag. 38
SL	ODPORN PLOŠČA IZ KALJENEGA STEKLA	pag. 42

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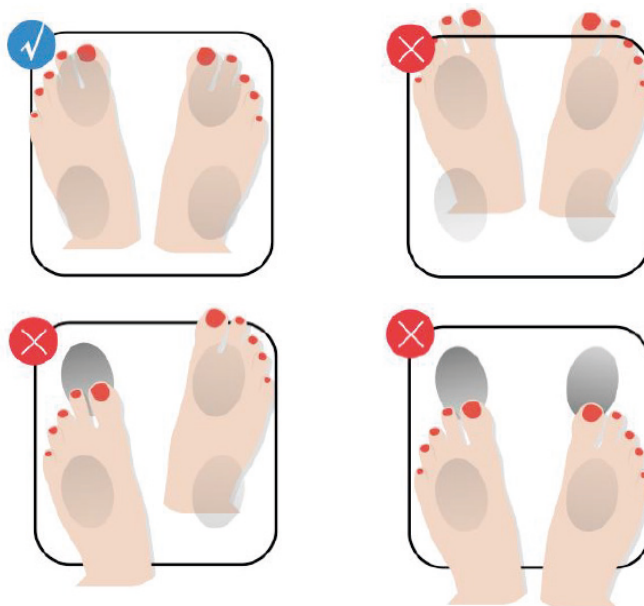**assistenza@beper.com**

## **General safety advice**

### **Read the operating instructions carefully before first use of the appliance.**

Read the following instructions carefully, since they give you useful safety information about installation, use and maintenance and help you avoid mishaps and possible accidents.

Remove the packing and make sure that the appliance is intact, with special attention to the supply cable. The elements of the packing (plastic bags, polystyrene, etc) do not have to be accessible to the children, in order to preserve them from danger sources.

It is recommended not to throw these packing material in the domestic waste, but to deliver them to the appropriate station of collection or destruction refusals, asking eventually information to the manager of the sanitation department of your own municipality.

This appliance is to be destined to its special purpose only.

Any other use is to be considered as improper and thus dangerous. The manufacturer cannot be responsible for any damages deriving from improper, wrong or incautious use.

To avoid any heating risk, it is recommended to unroll completely the supply cord and to unplug the appliance, when not used.

Before cleaning, always unplug the appliance.

To protect against electric shock, do not immerse cord, plug or any parts of the appliance into water or other fluids.

Do not hold the appliance with wet hands or feet.

Do not hold the appliance with naked feet.

Do not expose the appliance to atmospheric agents such as rain, wind, snow, hail.

Do not use the appliance near explosives, highly flammable materials, gases, burning flames, heaters.

Do not allow children, disabled people or every people without experience or technical ability to use the appliance if they are not properly watched over. Close supervision is necessary when any appliance used by or near children.

Make sure that children do not play with this appliance.

If you decide not to use the appliance any longer, make it unserviceable by simply cutting the supply cable after disconnecting the plug from the socket.

**Also make any dangerous parts unserviceable, especially for children who could play with the appliance.**

## Warnings

Do not place the body scale on slippery surfaces.

Make sure to place both feet at the center of the glass panel during weighting.

The panel of the body scale is made of tempered glass, make sure to avoid violent shocks and/or accidental drops, as the broken glass can cause wounds.

Store the body scale in a dry environment.

Do not place the body scale in water.

Use a soft cloth or a mild detergent to clean the outer surface of the scale. Do not use organic solvents or other chemicals.

When batteries have run out, replace them with same type, new batteries.

When not in use for long periods of time, remove batteries to avoid accidental leakage of corrosive liquid.

During weighting, stand barefoot and facing the display.

Set all the correct parameters for genre, height, weight and other indicators, to obtain accurate weight results.

### Check settings on your device

Before use, make sure that the device that will be paired with the body scale has the following specifications:

Apple devices: Bluetooth 4.0 or above; iOS 9.0 or above

Android devices: Bluetooth 4.0; Android 4.3 or above

### Placing and replacing of batteries

Put the scale upside down. Open the battery case found at the base of the scale.

Insert the batteries, making sure to respect the indicated polarity.

Close the battery case.

### Operation

Once the batteries have been installed, the body scale is ready to be used.

In case it is not paired with the App by Bluetooth connection, the item works as a regular body scale, and it can be used to assess body weight.

**Important: Beper declines all responsibility in case of malfunction and/or failed updating of the App, because this App is property of VTrump Tech Co.Ltd.**

**For further information on the App, visit the dedicated section in the App Store/Play Store.**

**Note: Beper has no control on sensitive data collected by the App, as these are entirely managed by third parties, and it has no access to the data collected through the use of the App.**

**Read the privacy statement related to the App found in the App Store/Play Store.**

### Downloading and setting of the App

The App can be downloaded through the QR code provided here on the side

It also possible to download the App following the steps here below:

Activate Bluetooth on your smart device

Make sure to connect your device to a wi-fi or 3G/4G network

Open the App store on your device and look for "ScaleUp"; proceed with download.

Once the App has been downloaded it will be possible to register the users and customize profiles based on personal information.

ScaleUp is able to assess 21 different body indicators: body shape, nutrition, weight, BMI, Body fat %, Muscles weight, Visceral fat, Water level, Bone structure, Basal metabolic rate, Body fat weight, Bone muscles, Water content weight, Protein content weight, Ideal weight, Weight control, Body fat control, Muscular mass control, Body age, Protein content %, Lean body mass.

#### **Instructions for use**

After setting your personal data and information, keeping the App running on your device step on the body scale barefoot, making sure your feet are correctly positioned on the glass panel Fig. A.

Wait until the following data are displayed on your device:

#### **Body fat content**

Related to the content of body fat calculated over body total weight; obesity should be assessed through body fat content rather than body weight.

#### **Basal metabolic rate**

Basal metabolic rate (BMR) is the quantity of energy that is used throughout the day to keep vital functions running (breathing, heartbeat, heat conservation, etc.).

#### **Water content**

Related to the body water content; with average body fat, body water content should be as follows:

Women 45% - 60%; Men: 55% - 65%

#### **Protein**

Main substances that make up human organs and tissues, their role is very important in all human life activities. Bone mass

In this case, bone content is related to the bone salt content, which is the mineral content of the bones. Bones are constantly undergoing an active metabolism that provides the calcium the body needs, always creating new bone mass. Being too thin or not getting enough exercise is closely related to bone loss.

#### **Visceral fat**

This is the fat that accumulates around intestines, it is easily accumulated and decomposed. Visceral fat is considered the most important, as well as one of the disease risk assessment factors.

Please note that the assessment of all these indicators can easily vary due to many different factors, it is thus advisable to assess weight always at the same time during the day, preferably before exercising and away from meals.

#### **Cleaning and maintenance**

The body scale can be cleaned with a soft cloth slightly dampened with warm water.

Do not use thinners, or other abrasive detergents, as these might damage the surface of the body scale.

Do not immerse the scale in water.

Protect the scale from accidental infiltration of liquids.

Keep hair lotions, make up or other soaps and lotions away from the body scale, as these can damage the appearance and/or the inner circuit of the scale.

#### **Tecnical data**

Power source: 2xAAA batteries (not included)

**For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.**

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment.

The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

### **GUARANTEE CERTIFICATE**

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

**The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance**

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

### **GUARANTEE CONDITIONS**

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product

Contact your distributor in your country or after sales department Beper.

E-mail [assistenza@beper.com](mailto:assistenza@beper.com) which will forward your enquiries to your distributor.

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## YOU MIGHT ALSO LIKE

### Rechargeable Percussion Massager

Cod.: P302MAS001

Ideal for the deep and targeted treatment of muscular soreness and pain.

6 adjustable intensity levels

4 interchangeable massaging accessories

Long lasting lithium battery

Ergonomic and anti-slip handle

Low noise and light weight for easy use

## YOU MIGHT ALSO LIKE

### Vacuum Blender

Cod.: BP.620

Reduces oxidation and keeps the smoothie fresh and nutrient

Accessories: 800 ml BPA free travel container,

lid with date indicator , vacuum tube

2 preselected modes: Vacuum Smoothie and Smoothie

Adjustable speed and Pulse function

Touch screen control panel

Steel blades

1L BPA Free Jar

## **BEPER SRL**

Via Salieri, 30 - 37050 - Vallese di Oppeano - Verona  
Tel. 045/7134674 – Fax 045/6984019  
[beper.com](http://beper.com)

## **CUSTOMER CARE BEPER**