beper

- SPAZZOLINO RICARICABILE SONIC
MANUALE DI ISTRUZIONI
- RECHARGEABLE SONIC TOOTHBRUSH
USE INSTRUCTIONS
- BROSSE À DENTS RECHARGEABLE SONIC
MANUEL D'INSTRUCTIONS
- WIEDERAUFLADBARE SONICZAHNBÜRSTE
BETRIEBSANLEITUNG
- CEPILLO DENTAL RECARGABLE SÓNICO
MANUAL DE INSTRUCCIONES

Cod.: 40.913

ITALIANO	pag. 3
ENGLISH	pag. 8
FRANÇAIS	pag. 12
DEUTSCH	pag. 17
ESPAÑOL	pag. 21

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e-mail: assistenza@beper.com

Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.

Warning

Check if the voltage indicated on the bottom of the charger correspond to the local mains voltage before you connect the appliance.

The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.

Always have the charger replaced with one of the original type in order to avoid a hazard.

Do not use the charger outdoors or near heated surfaces.

If the appliance is damaged in any way (brush head, toothbrush handle, charger), stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Customer Care Centre in your country (see chapter 'Guarantee & service').

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.

Caution

Do not clean the brush head, the handle, the charger, the charger cover in the dishwasher. If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush. Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use.

The Sonic toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.

If you have medical concerns, consult your doctor before you use the Sonic.

This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.

The Sonic toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.

Do not use other brush heads than the ones recommended by the manufacturer. If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.

The Sonic

- A Hygienic travel cap
- B Brush head
- C Interchangeable brush head
- D Handle with soft grip
- E Power on/off button and mode selector
- F Brushing modes indicator
- G Charger indicator
- H Charger base

Attaching the brush head

- 1 Align the brush head so the bristles face the front of the handle.
- 2 Firmly press the brush head down on the metal shaft until it stops.

Using the Sonic

Brushing instructions

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle towards the gumline.
- 3 Press the power on/off button to switch on the Sonic.
- 4 Apply light pressure to maximise the Sonic's effectiveness and let the Sonic toothbrush do the brushing for you.
- 5 Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth.

Continue this motion throughout your brushing cycle.

- 6 Fig. 2 Begin brushing in section 1 (outside of upper teeth) and brush for 30 seconds before moving to section 2 (inside of upper teeth). Continue brushing in section 3 (outside of lower teeth) and brush for 30 seconds before moving to section 4 (inside of lower teeth).

Sonic make a pause every 30 seconds to moving to next section.

- 7 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush switched on or off, as you prefer.

After 2 minutes Sonic will auto shut off.

Brushing modes

Clean mode


Standard mode for superior teeth cleaning.

Sensitive mode

Gentle, yet thorough cleaning for sensitive gums and teeth.

Massage mode

Gentle gum stimulation.

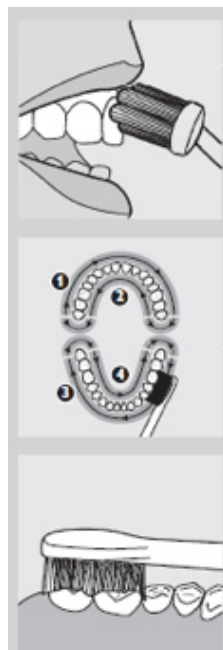Cleaning

Do not clean brush heads, the handle, the travel charger, the charger cover and in the dishwasher.

Toothbrush handle

1 Remove the brush head and rinse the metal shaft with warm water.

Do not push on the rubber seal round the metal shaft with sharp objects, as this may cause damage.

2 Use a damp cloth to wipe the entire surface of the handle.

Brush head

1 Rinse the brush head and bristles after each use.

2 Remove the brush head from the handle and rinse the brush head connection at least once a week with warm water.

Storage

If you are not going to use the appliance for a long time, remove the mains plug from the wall socket. Then clean the appliance and store it in a cool and dry place away from direct sunlight. You can use the cord wrap to store the mains cord neatly.

TECHNICAL DATA

Power 2W

Power supply 100-240V ~ 50/60Hz

For any improvement reasons, Beper reserves the right to modify or improve the product without any notice

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment. The crossed out “wheeled bin” symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

Guarantee conditions

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- Repairs are professionally carried out by qualified persons who are appointed by the supplier.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) is excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

Technical assistance

Even after the guarantee period we will always pay attention to repairs for defected appliances. For technical assistance and/or repairs after the guarantee period you can directly contact the below address:

CONTACT YOUR DISTRIBUTOR IN YOUR COUNTRY OR AFTER SALES DEPARTMENT
BEPER. E-MAIL assistenza@beper.com
WHICH WILL FORWARD YOUR ENQUIRIES TO YOUR DISTRIBUTOR.

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